

Cat. No. BYD009-1ML

100 bp YEA Ladder DNA Marker III

Ready-to-use DNA ladder

Recommended storage condition: -20°C

| | |
|-----------------|--------------|
| Loading volume: | 5-10 µl/well |
| Size range: | 100-3000 bp |
| Concentration: | 120 ng/µl |
| Package size: | 120 µg/ml |

100 bp DNA Ladder Plus 1500 bp, 2000 bp and 3000 bp, ready-to-use

The Yeastern Biotech's 100 bp YEA Ladder DNA Marker III (BYD009-1ML) contains 13 discrete fragments ranging from 100 bp to 3000 bp. It is supplied in a ready-to-use format containing blue tracking dyes, and mixed at varied concentrations.

For best results, load 5 µl of the marker per well. It allows sizing and concentration estimation of DNA fragments on agarose gels generated by PCR or restriction digestion.



www.yeastern.com


yeastern@yeastern.com

t / 886-2-26953922 f / 886-2-26953979

Cat. No. BYD009-1ML

100 bp YEA Ladder DNA Marker III

Ready-to-use DNA ladder

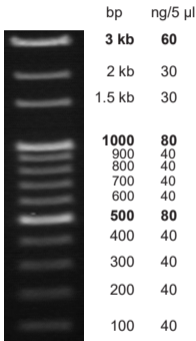

Loading volume: 5 μ l/well
1.2%, 1 \times TAE agarose gel



Yeastern Biotech Co., Ltd.